 **SMART SCALE**  
BILANCIA SMART IMPEDENZIOMETRICA PESAPERSONE

MANUALE D'USO  
USER MANUAL  
BEDIENUNGSANLEITUNG  
MANUEL D'UTILISATION  
MANUAL DE USO

# INDEX

|    |       |          |
|----|-------|----------|
| IT | ..... | P. 03-07 |
| EN | ..... | P. 08-12 |
| DE | ..... | P. 13-17 |
| FR | ..... | P. 18-22 |
| ES | ..... | P. 23-27 |

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
## 01 INSERT THE BATTERIES

1. Open the battery compartment on the back of the scale;
2. Insert the batteries following the polarity as indicated;
3. Close the compartment.

### Important:



- Always use the correct size and type of batteries.
- Do not mix new and used batteries or batteries of different types or brands.
- Before installation, clean the battery contacts and those of the scale.
- If 'Lo' appears on the display, be sure to replace all batteries in a set at the same time.

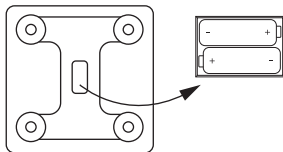The batteries required for the scale to operate are 2 AAA batteries.

## 02 SCALE POSITIONING

- Use the scales on a flat, stable surface.
- Avoid uneven surfaces such as carpets or floors with uneven tiles.
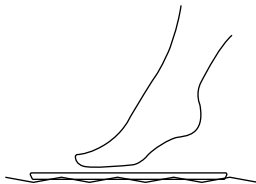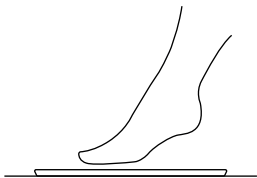
### 03 DOWNLOAD THE APP


1. Turn on Bluetooth;
2. Ensure that your phone/tablet is connected to the internet.

Use your device to scan the QR Code to download the OKOK International app, or search for the app on the App Store or Google Play.

### 04 REGISTER AND COMPLETE WITH PERSONAL INFORMATION



1. Register by mail;
2. Select the unit of weight and enter your personal information.
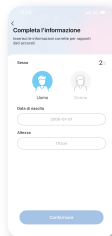
## 05 SYNCHRONISE THE DEVICE AND MEASURE THE WEIGHT

1. Click on the 'Add device' button or press the '+' button to select the scale;


2. Step onto the scales with bare, dry feet and wait for matching.
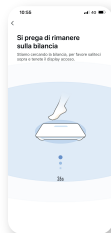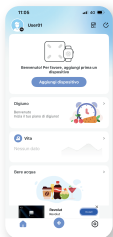
## 06 TRACK PROGRESS AND SET GOALS

1. Go to the "Weight Trend" page;

2. View progress and trends.

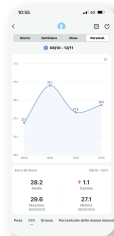- Display: All data displayed when switched on
- Capacity: 3 kg - 180 kg
- Division: 1D = 100 g / 0,2 lb
- Units: KG / LB / ST (Kilograms / Pounds / Stone)
- Voltage: 4,5 V
- Battery: 2 x AAA batteries
- Low battery indication: 2,4 V  $\pm$  0,2 V, displays 'Lo' and switches off
- Current: < 25 mA
- Standby current:  $\leq$  10  $\mu$ A
- Overload indication: overload exceeding 2,5kg, for a total of 182,5kg displays 'Err'
- Best operating temperature: 5 ~ 35 °C
- Storage temperature: -10 ~ 60 °C
- Humid environment:  $\leq$  90% RH

## Instructions for the OKOK International App

### Instant weighing

Open the OKOK APP, it will remind you to turn on Bluetooth. Turn on Bluetooth and the automatic pairing process will be completed automatically. When the user starts weighing themselves, the weight data will be synchronised with the APP. If Bluetooth is disconnected, the weight data will be saved in the scale's memory and then synchronised with the APP once Bluetooth is reconnected.

### Manual entry of weight data in the APP

The user can manually enter weight data in the APP to track progress.

### Current Weight

The Current Weight interface displays information about the user's and family members' weight, and generates a report providing health advice based on the measurement. It includes the Current Weight interface, Weight List interface, Report by Date, and Notifications.

### Weight trend

Includes weekly, monthly and quarterly reports showing weight changes using curve graphs.

### Notice

For the body fat scale, there are 8 data points: weight, body fat, visceral fat, muscle, bone, BMR (Basal Metabolic Rate), age, BMI (Body Mass Index).

Scale management is divided into two main areas: scale management and connection management.

- Scale functionality:

This includes managing weighing operations and the process of pairing, associating and disassociating devices.

- Connection management:

This includes checking the operating status and managing the Bluetooth connection.

### **Account and user**

Includes account management and user management:

- Account management:

Includes registration, login, password change, password recovery, links to Facebook, Twitter, etc.

- User management:

Allows management of user and family member information (nickname, birthday, height, gender, etc.) all associated with the same registration account.

### **Warnings**

1. Keep the scale away from wet floors. Do not use the scale if the glass surface is wet to avoid slipping.
2. Step onto the center of the scale.
3. The scale platform is made of glass, avoid sharp or rough objects that could scratch the surface. Excessive forces, sudden shaking, or unexpected falls can break the glass and injure the user.
4. The Bluetooth scale is not waterproof. Keep the surrounding environment dry and avoid submerging it in water.
5. Use a soft cloth with neutral detergents to clean the surface, avoiding contact with organic solvents or chemicals.
6. When the batteries are dead, remove them to prevent any leaking liquids from corroding the scale.

### **Notices**

1. Step on the scale correctly.
2. Stay on the scale for no less than 10 seconds if measuring body fat.
3. Make sure your gender, height, and weight are correct.

### **Warranty**

The digital weighing scale is covered by a guarantee in accordance with the law.

The warranty does not cover damage due to wear, accidental breakage, incorrect use of the appliance due to non-compliance with the user manual or negligence in use and maintenance.

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